
QUICK START

1. Remove packaging and place your SunTouch Plus on a level surface.


2. Plug unit into an electrical socket.

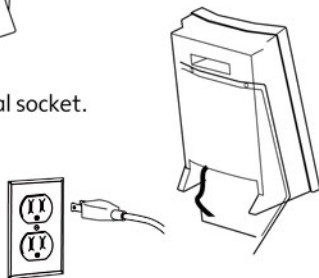3. Using the control panel adjust unit to desired setting.
(See right for control instructions)

4. Place unit 12 to 16 inches away and direct light toward face for 20 to 30 minutes per day.

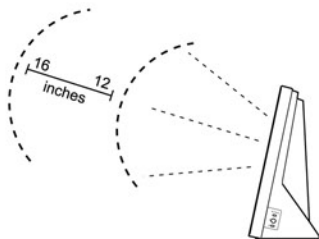
CONTROL PANEL

Turn the unit on by pressing the **on/off** button. The indicator will light up.

Use the **timer** button to select your desired duration of use.


To turn on the negative ion generation, you may press the **ion** button. The indicator will light up when ion option is working.

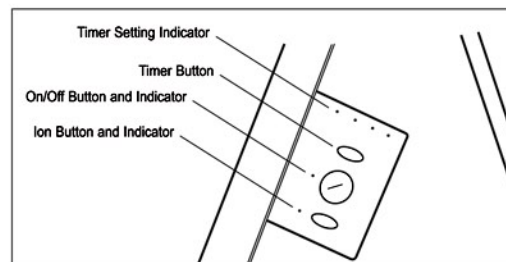REPLACEMENT PARTS

Replacement parts for Sun Touch Plus:
SunTouch Plus Bulb,
Part Number: F4040B

Contact your NatureBright Dealer
for ordering information.

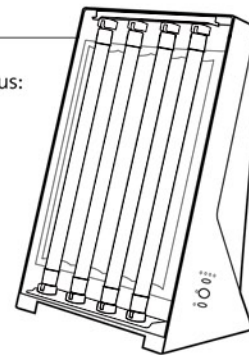NatureBright®


SunTouch

PLUS

LIGHT AND ION THERAPY

User's Guide

- Setup
- Operation
- Specifications
- Warranty

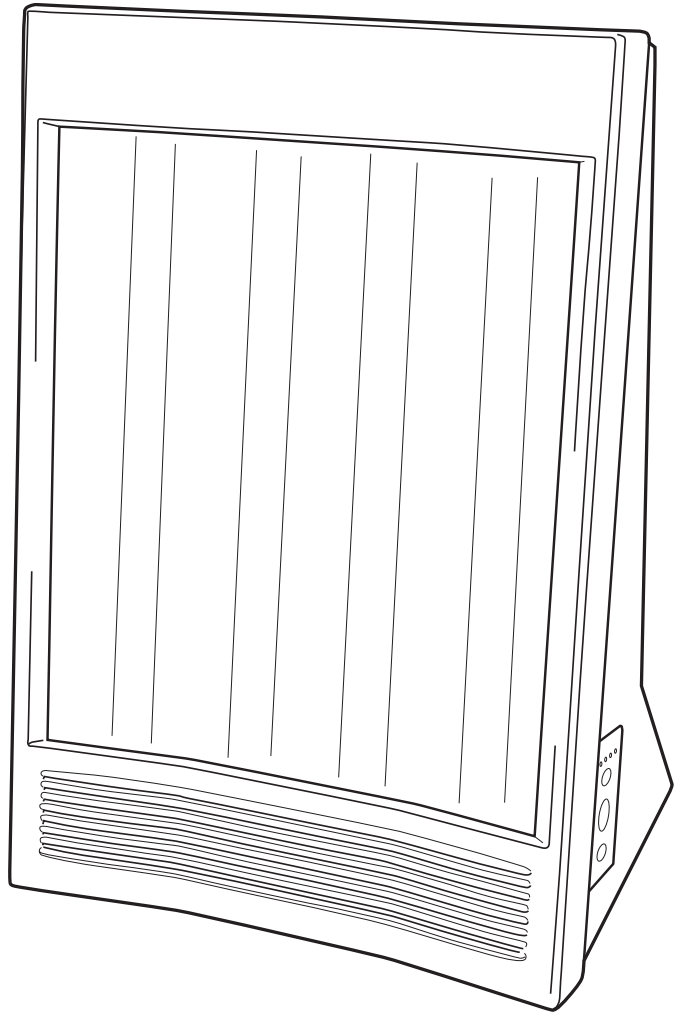CAUTION: Read manual carefully before operating.

Nature Bright Company 2600 MICHELSON DRIVE, 17TH FLOOR, IRVINE, CA 92612 USA

© Nature Bright Company 2007

Contents

SunTouch Plus is powered by SkyEffect™ Technology.....	3
SunTouch Plus Negative Ion Therapy.....	3
Product Specifications.....	4
Unit Inspection.....	4
Placing the Unit.....	5
Operating Instructions.....	5
Post Purchase Service.....	6
Frequent Asked Questions.....	7
Replacing the Bulb(s).....	8

WARNING – HIGH VOLTAGE

CAUTION!

Rise of Fire-Use only with 8 Watt Bulbs

CAUTION! – This equipment should be inspected frequently. Dirt that collects on the unit should be removed regularly to prevent fire.

CAUTION! – For AC outlet installation only. Never connect to an extension cord or other power source. Improper installation could result in electrical shock or short.

SunTouch Plus is powered by SKY EFFECT™ Technology

NatureBright developed Sky Effect, a new light therapy modality, to leverage a recent breakthrough discovery of a third light receptor in the human eye. This retinal receptor is sensitive to light frequencies. It monitors or “reads” wavelengths, particularly of blue light, sending signals to the brain’s hypothalamus, the area of the brain where sleep, hormonal, and emotive functions are regulated. As this receptor sends signals to the hypothalamus, it regulates our body clock, that cellular sensitivity that keeps us in tune with the earth’s natural time and seasonal cycles.

The more in sync our body clock is with optimal natural cycles, called the Circadian rhythm, the better it will work to stimulate our overall well-being. We’ll sleep well, want to eat at regular times, and experience emotional stability. If our clock is out of sync, which can happen from illness, night shift work, jet lag, or winter light deprivation, we can experience disturbing symptoms such as insomnia, and other problems that can range from daytime drowsiness to depression.

Sky Effect Technology, used in the SunTouch Plus and other NatureBright products, is entirely safe, gentle, and soothing. It can be used, according to directions, without any negative side effects to produce a sense of long term health and emotional and physical well being.

SunTouch Plus Negative Ion Therapy

Odorless, tasteless, microscopic negative ions are found in the air, naturally. Research shows negatively charged ions can lift your spirits. A recent scientific study on depression, mood, and Seasonal Affective Disorder (SAD) found that 58% of patients studied experienced a significant improvement when they were exposed to streams of high-density ions.

Flowing negative ions into the air also helps asthmatics and those with allergies. It eliminates tiny irritative particles such as pollen, mold spores, and other irritants from the air.

Safe & Easy to Use in the Home or Office

Negative ions are safe to use in the home or office. As easy to operate as an air purifier or home humidifier, simply running the SunTouch Plus can lift your spirits, mental alertness, energy levels, and mood—safely and gently with no negative side effects.

Product Specifications

Electrical:

- 120V, 0.4 A, 60Hz
- 45 Watt power consumption
- 4X8 Watt T5 Tubes

Output:

- 10,000 LUX

Size:

- Height: 13.375" (34 cm)
- Width: 8.125" (2.6 cm)
- Depth: 7" (17.8 cm)
- Weight: 7 lbs (3.2 kg)

Warranty:

- Length of warranty varies depending on model.
- See warranty card for specific details.

Unit Inspection

Read all safety guidelines before using your SunTouch Plus Light.

These safety instructions must be observed in order to avoid potential hazards that could result in personal injury or damage.

Never use damaged or worn cords or plugs. The use of damaged or worn power cords or plugs could result in electric shock, burns, or fire. Do not place heavy objects on top of the cord. When unplugging the unit, grasp the plug directly to avoid damage to the cord. Never pull on the power cord in an attempt to remove a plug from the socket.

Unplug the unit before servicing to eliminate the risk of electrical shock.

SunTouch Plus is designed to allow the front cover to be removed and the bulbs to be changed safely.

Any other disassembly will void the warranty.

The lamps get extremely hot. Allow five minutes for the lamps to cool before handling.

To avoid damage to the light box, do not allow any liquid contact any part of the light.

Placing the Unit

Place the light near to you at eye level.

For the best results, place your SunTouch Plus so its light bathes your face. Many people use SunTouch Plus on a counter, table or exercise while using it. You can read, eat or work at the computer, watch TV or even when putting on makeup. Use in a well-lit room to avoid eyestrain.

Turn towards the light.

The proper hormonal response to light therapy is through the retina of the eyes, not through the skin. However, you do not need to stare at the light. If you are in front of the light and it is close enough, only an occasional glance is necessary.

This is not a tanning light. SunTouch Plus creates specific wavelengths of light and reduces UV exposure. The light is flicker-free to reduce eyestrain.

Operating Instructions

Powering your SunTouch Plus

The SunTouch Plus has been designed with a six and half foot cord that gives you freedom to move the unit to the most preferred position.


After plugging the cord into a standard electrical outlet, turn the unit on. The SunTouch Plus is turned on/off by a switch located on the right hand side of the unit. The SunTouch Plus comes with two switches. One controls the lamp, the other the ionizer. Each switch works independently giving you more options. Once you are finished with your session, you can turn off the light and still enjoy the benefits of the ionizer throughout the day.

Results

Most people feel a difference right away and may obtain measurable results in as little as four days.

Optimizing your light Therapy

You have to sit fairly close to have an effective session. For optimal results sit twelve (12) to sixteen (16) inches away. If you have the SunTouch Plus with the added benefit of ionization, you can sit up to eighteen inches away and achieve similar results. For best results, allow the light to bathe your face.

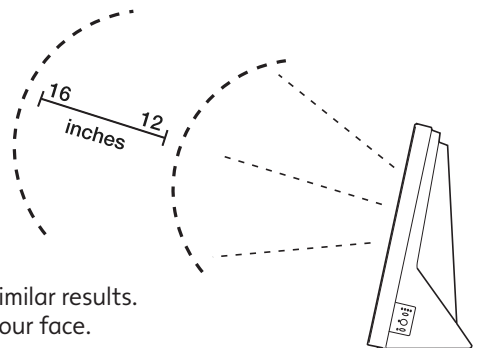Parameters for Using the SunTouch Plus®

If you encounter following problems, related to Circadian Rhythm disturbances:	Use these guidelines for Morning Light Therapy:	Use these guidelines for Evening Light Therapy:
Winter SAD	30 min upon awakening	
Low Mood	30 min upon awakening	
ASPS: falling asleep too early		30 min between 8 pm to 9 pm
DSPS: falling asleep too late	30 min between 7 am - 8 am	
Irregular sleep-wake phase	30 min between 6 am - 8 am	
Shift Work Fatigue		30 min at home before the start of your shift
Jet Lag	Two 30 min morning sessions on consecutive days after eastward travel	Two 30 min evening sessions on consecutive days after westward travel

Session time – First use

Leading researchers recommend thirty minutes of 10,000 LUX, light therapy powered treatment. All SunTouch Plus powered lights provide 10,000 LUX, but at varying distances.

Reducing session time

After a week or two of using the SunTouch Plus, once your desired results are reached, you can reduce your session time from 30 minutes to 15 or 20 minutes per session.

Replacing Bulbs

The bulbs inside the light box usually do not burn out, but rather decrease in intensity overtime. We recommend that you replace your lamps every two years. Replacement bulbs can be ordered online or by calling NatureBright's Service Department.

Post Purchase Service

In the event that your SunTouch Plus does not function properly or one of the bulbs stops working, please call NatureBright's Service Department at 1-800-622-0231.

Frequent Asked Questions

1. Do I have to take off my glasses when I use the unit?

Not usually, but if you have transitional lenses that darken when you go outside, you will have to remove them during your session.

2. Do I need to clean the ionizer?

The ionizer doesn't require internal cleaning. With repeated use you may get some build up around the vent area. This is easily wiped away with a dust cloth.

3. Are there any major side effects?

No. The SunTouch Plus is a safe product.

4. Can I get sunburn from using SunTouch Plus?

No. While typical fluorescent bulbs produce UV light, SunTouch Plus eliminates almost all measurable UV light.

5. Shall I keep my eyes open for optimal results?

The SunTouch Plus light works through the retina of the eyes, not the skin. This is not a tanning light. In fact, SunTouch Plus creates specific wavelengths, which produce effective results while eliminating UV hazards and reducing eyestrain. Keep in mind that you don't have to stare at the light. As long as you are close to the light and glance at it on occasion will suffice.

6. What are the suggested duration times?



Duration times vary according to your needs. If your physician has recommended a plan for you, follow that plan. Most people will obtain best results using the light in the morning. Sessions can be as short as fifteen minutes. A majority of people achieve good results with 20 to 30 minutes of usage. You may want to experiment to find your own personal optimal duration.

7. Should I use SunTouch Plus every day?

You should use it daily, especially during the winter months. Many people find that they don't need the light during the summer. Those working long hours indoors may find the SunTouch Plus to be effective year round. About 50% of the time, those who need to use morning light find that they also benefit from evening light. If this describes you, make sure to sit about 12 inches farther away from your light in the evening to help prevent conflicting signals with morning light.



Replacing Bulbs

1. Lay the SunTouch Plus on it's face and remove the four (4) face-plate screws.

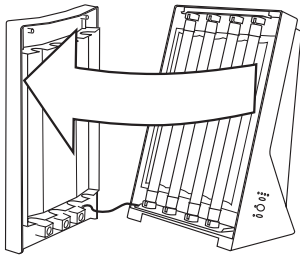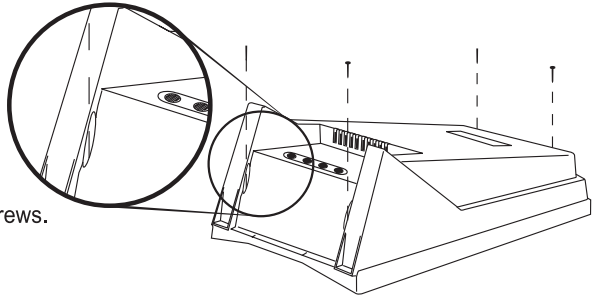
2. Remove cover being careful not to pull on the ground wire attached to faceplate. Lay faceplate aside.

3. Remove bulbs by twisting to line-up pins with the opening (1). Start by tipping out top of bulb, then pull bottom out (2).

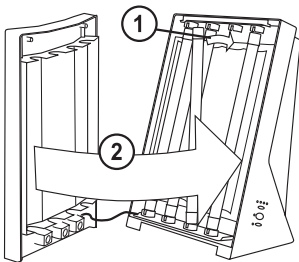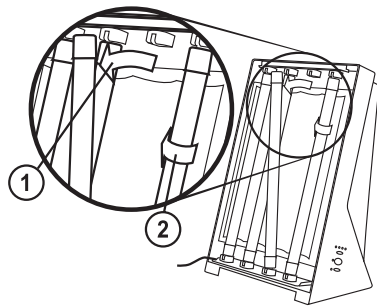4. Insert new bulb into sockets. Twist 1/4 turn to lock. Replace cover.

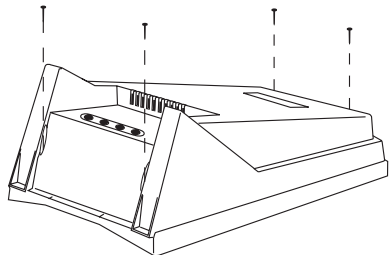5. Lay the unit on face and replace four (4) screws. Tighten to snug fit.

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